2017-01
801514
1.0
M

 TRACOE medical GmbH
Reichelsheimer Straße 1 / 3, 55268 Nieder-Olm / Germany
Tel.: +49 6136 9169-0, Fax: +49 6136 9169-200
info@tracoe.com, www.tracoe.com

TRACOE®
twist

CE
0 1 2 3

**TRACOE® *twist* Tracheostomy Tubes
with Low Pressure Cuff and Atraumatic Insertion**

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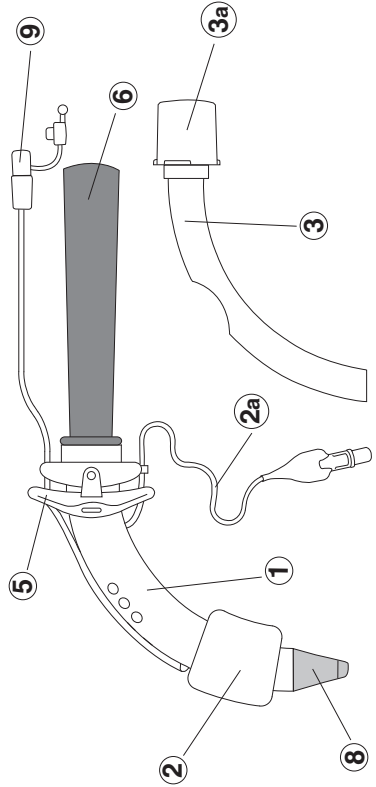
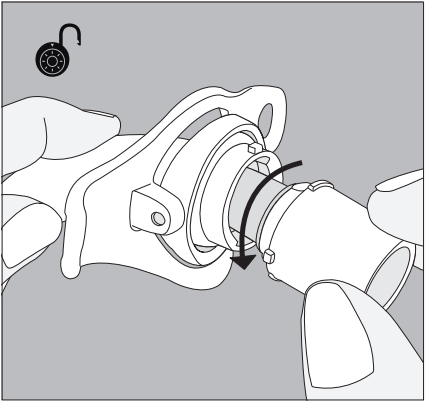
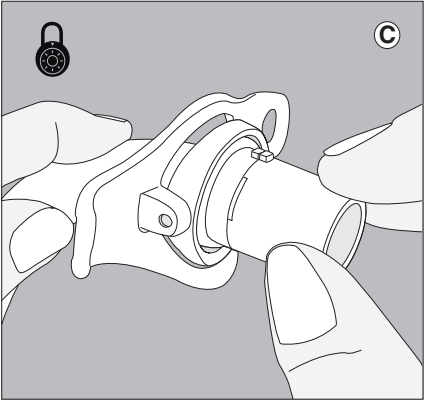
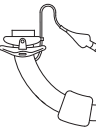
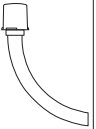
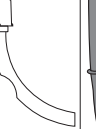
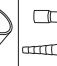
Instructions for Use for TRACOE® *twist* Tracheostomy Tubes with Cuff and Atraumatic Inserter

NOTE:

- Please read the instructions for use carefully. They are part of the described product and must be available at all times. For your own and your patients' safety, please observe the following safety information.

- The illustrations to which the text refers can be found on the (foldout) illustrated pages at the beginning of these instructions. The numbers and letters in brackets refer to the respective illustrations and product components of the tracheostomy tube. The symbols and icons used are explained on pages 1 to 5.

1. General Information

CAUTION: The product may only be used by **trained** doctors and nursing staff. The attending doctor must determine the correct size and type.

Intended Purpose: This medical device is a tracheostomy tube with an inserter for a PDT (percutaneous dilation tracheostomy) using a Seldinger wire $\varnothing 1.27 \pm 0.04$ mm (e.g. REF 518, not included in the delivery) and reinsertion.

Functional Description: The tracheostomy tube ensures that respired air passes through the tracheostoma.

- The inflated cuff separates the upper from the lower respiratory tract, preventing the flow of air from the lungs into the mouth and nose and vice versa. This ensures that the patient can breathe **only** through the tube. The cuff is inflated with air using the inflation line. When inflated, the cuff seals the space between the trachea and the outer cannula wall. This allows efficient ventilation with ventilation equipment and prevents subglottic secretions from entering the lower respiratory tract.
- Fenestrated cannulas allow some of the air needed for speaking to reach the upper respiratory tract through the fenestration. Airway resistance in the upper respiratory tract is also reduced.
- The atraumatic inserter bridges the gap in diameter.

- See below for a more detailed functional description.

Maximum Useful Life of the Tracheostomy Tube:

- With cuff – 29 days from first use (see chapter 5 "General Precautions"). This maximum period of use also includes any periods when the tracheostomy tubes or inner cannulas are not in use.

Product for Single-Patient Use: Multiple use by one and the same patient during the period of use is permitted.

Single-Use Product: The atraumatic inserter (6, 8) is a single-use product. Do not clean or reuse.

Product Info Card: This product is supplied with a product info card with two detachable labels. The labels show product details. Keep the info card separately, to facilitate reordering. The detachable label can be affixed to the patient's medical records, for example.

2. General Description

The product consists of an outer cannula (1) with a cuff (2) or without a cuff, and an exchangeable inner cannula (3) with a 15 mm connector (3a) or grooved low profile connector. The outer and inner cannulas are connected and disconnected with a twisting motion (C).

The cuff is inflated or deflated using the inflation line (2a).

The outer cannula is made of radiopaque polyurethane and is attached to a neck flange (5) that can swivel on 2 axes (cardanic suspension).

The inserter is already mounted in the outer cannula and has a foldable silicone sleeve at its tip (8).

3. Indications

Due to the pre-installed atraumatic inserter, the tubes can be used for percutaneous dilation tracheostomies in combination with the TRACOE experc dilation set REF 520 or for tube change. The Seldinger technique must be applied in both cases.

The products are suitable for patients where access to the respiratory tract is required via a tracheostoma with the trachea sealed. The fenestrated tubes with cuff facilitate speech in patients with an intact larynx.

The tubes with a suction line (9) are for use in patients where suction of subglottic secretion is indicated.

4. Contraindications

4.1 Absolute Contraindications

- Not suitable for emergency tracheostomy/cricotomy.
- Existing infections in the tracheostomy area,
- Existing malignancy in the tracheostomy area,
- Unstable fracture of the cervical spine,
- Unclear identification of anatomical marker points, and
- unsuitable anatomical conditions.
- Do not use the occlusion cap (**10**)/speaking valve in laryngectomised patients (patients without a larynx) – risk of suffocation!

4.2 Relative Contraindications

Relative contraindications are cases where the risks and benefits of the procedure must be carefully evaluated:

- An enlarged thyroid,
- Previous surgical intervention in the neck and throat area (e.g. thyroidectomy),
- Strong tendency to bleeding, e.g. during treatment with anticoagulants,
- Unusually deep trachea, e.g. in obese patients, where an extra-long tube may be required,
- Paediatric use.

5. General Precautions

- For the first percutaneous dilation tracheostomy, always keep an identical set at hand as a replacement.
- You are strongly advised to keep a ready-to-use spare tube and several spare inner cannulas at the patient's bedside. Store in a clean and dry state.
 - Check the tube for integrity and correct function prior to use/insertion. Check, for example, for the absence of obstruction, cuff seal, correct and stable fit of the inner cannula inside the outer cannula, absence of kinks, stable connection between the tube and the neck flange, etc. The cuff material must not be brittle. If the product is damaged, replace with a new product.
 - Do not use force on the tracheostomy tube, as this may damage or break it. If the connections at the 15 mm connector are tightly fixed, always use a disconnect wedge approved for tracheostomy tubes.
 - Keep the 15 mm connector (**3a**) clean and dry.
 - When changing the inner cannula, always ensure that the inflation line (**2a**) of the cuff is not positioned between the inner and outer cannulas, as it may get trapped and damaged.
 - When repositioning the patient in bed, ensure that the patient does not lie on the pilot balloon (**2a**), as this would greatly

increase the cuff pressure and may damage the trachea.

- During mechanical ventilation and frequent changes of the patient's position or manipulation of the tube, the inner cannula may become separated from the outer cannula.
- With the use of high ventilation pressures, leakage between the inner and outer cannulas may occur in isolated cases.
- The use of tracheostomy tubes can lead to pressure points, necrosis of the skin on the neck and skin irritation (e.g. moisture). To avoid this, we recommend placing a pad underneath the neck flange.
- To avoid damage to the material, the cuff should not be in contact with lidocaine-containing aerosols or any ointments.
- All components of the cuff inflation system must be kink-free and not taut when testing the cuff pressure. The hand-held manometer may otherwise give a false pressure reading.
- Thin cuffs are to some extent permeable to water vapour. Therefore, condensate may accumulate in the cuff. Where quantities are small, this is of no significance. However, if larger quantities of condensate are inadvertently aspirated into the inflation line, it is no longer possible to measure and adjust the cuff pressure correctly and the tube needs to be replaced. Before removing the tube, remove as much air or water from the cuff as possible with a syringe.
- When using together with other medical devices, observe the respective instructions for use. If in doubt, contact the manufacturer.

6. Warnings

- Do not use if the sterile packaging is not intact.
- The atraumatic inserter supplied with the tracheostomy tube must be used only for the corresponding tube of the appropriate size. The size is clearly marked on the inserter.
- Never use fenestrated inner cannulas for ventilation.
- Use the occlusion cap only with fenestrated tubes (outer + inner cannula) and an unblocked cuff.
- Do not use a speaking valve during sleep.
- When using laser or electrosurgical treatment equipment, ensure that the tracheostomy tube is at a safe distance. There is a risk of fire and toxic gases and the tube may get damaged.
- Select the correct tube to ensure that the fenestration (if applicable) is positioned at sufficient distance from the stoma canal. Failure to do so exposes ventilated patients to the risk of emphysema. Granulation tissue may also form or airway resistance may increase if speaking valves or occlusion caps are used.
- If a fenestrated tube is used, there is a risk of increased formation of granulation tissue.

- The cuff pressure can change, for example due to altitude changes (e.g. in an aircraft), or if laughing gas is used as an anaesthetic or during the connection and disconnection of a hand-held manometer.
- Excessive cuff pressure poses a risk of permanent damage to the trachea.
- Excessive cuff pressure may give rise to cuff hernias.
- Insufficient cuff pressure poses a risk of aspiration.
- Ensure that the correct Luer connectors are used (e.g. **2a + 9**).
- During insertion and removal of the tube, irritation, need to cough or bleeding may occur.

7. Adverse Reactions

Pressure points, necrosis, skin irritation, granulation tissue, need to cough, difficulty swallowing, and bleeding.

8. Inserting the Tube

8.1 Preparing the Tube

For the first percutaneous dilation tracheostomy, always keep an identical *twist* set at hand as a replacement.

Carry out the following steps under sterile conditions.

1. Check the pack content for completeness (**D**).
2. Check the cuff (**2**) for leakage by inflating it. For this purpose, inflate the cuff with a hand-held manometer to a pressure of 50 cmH₂O (≈ 36.78 mmHg) and wait for 1 minute to see whether the cuff deflates. If the cuff is leak-tight, remove all air from the cuff with a syringe. Push the cuff upwards in the direction of the neck flange (**5**) to facilitate insertion.
3. Before using the tracheostomy tube, check that
 - the folding silicone sleeve at the distal end of the inserter fits seamlessly with the distal end of the tracheostomy tube
 - the guidewire passes freely through the guiding catheter.
4. To facilitate insertion, apply a pea-sized amount of the enclosed lubricating gel to the inserter and the tracheostomy tube.

8.2 Preparing the Patient

- Ensure that the patient is well pre-oxygenated before cannulation or recannulation.
- The patient should be supine with a cushion under the shoulders for the head and neck to be extended. The headrest of the patient's bed should be raised by 30 to 40 degrees.
 - Sedate the patient and monitor the vital functions.
 - Should complications arise during or recannulation, take appropriate safety measures (e.g. tracheal spreader, tube with

a smaller diameter) to allow the doctor to carry out short-term ventilation using translaryngeal intubation or a larynx mask.

8.3 Inserting the Tube after Dilation Tracheostomy

The following steps must be carried out by bronchoscopy:

1. Once the tracheostomised patient has been prepared for cannulation, insert the Seldinger wire and guiding catheter into the stoma canal. Now push the inserter over the unit (Seldinger wire/guiding catheter) as far as the safety stop.
2. Make sure that the proximal end of the guiding catheter always corresponds to the proximal marking of the Seldinger wire. In this way, the guidewire and guiding catheter will always precede the inserter.
3. Insert the unit of the recannulation tube or guiding catheter/Seldinger wire and the inserter with the tracheostomy tube into the trachea, and advance them until the neck flange touches the skin. We recommend holding the tracheostomy tube and the inserter together in one hand during insertion.
4. Retract the atraumatic inserter and the recannulation tube or guiding catheter together with the Seldinger wire, leaving the tracheostomy tube in place in the trachea. The easiest way of doing this is to hold the neck flange in place with one hand and remove the inserter by pushing against the neck flange with the thumb of the other hand.
5. Then insert the inner cannula. To secure the inner cannula, hold the neck flange with the fingertips and twist the connector of the inner cannula until it locks into place, i.e. until the markings are aligned with each other (**C**).

When inserting the inner cannula, ensure that the inflation line (**2a**) of the cuff is not positioned between the inner and outer cannulas, as it may get trapped and damaged.

Make sure that the silicone sleeve is still on the end of the inserter once it has been pulled out. If not, remove the radiopaque silicone sleeve from the tube or the respiratory tract.

6. For ventilation, connect the standardised 15 mm connector (**3**) of the inner cannula to a ventilation machine.
7. Only block the cuff (**2**) after ventilation via the tracheostomy tube has been established, otherwise there is a risk of asphyxiation or emphysema. Inflate the cuff using the pilot balloon (**2a**). Now check the cuff pressure to ensure that the cuff has not been damaged during insertion. Adjust the cuff pressure to the individual ventilation therapy and check at regular intervals. Typically, the pressure should be between 20 cmH₂O (\approx 15 mmHg) and 30 cmH₂O (\approx 22 mmHg).

Attaching a hand-held manometer causes a fall in cuff pressure.

This effect is more marked with small sizes.

8. Attach the neck strap (**12**) to the neck flange (**5**) to keep the tube in place.
9. Check the position of the tube (e.g. using bronchoscopy or diagnostic radiology, see chapter 5. "General Precautions") as well as its correct function (see chapter 6. "Warnings"). Carefully suction the trachea and tube to ensure that the airways are unobstructed.

8.4 Removing the Tube (see chapter 8.2)

Prepare to remove the tube as follows:

- Bend the patient's head slightly backwards and
- completely unblock the cuff of cuffed tubes (see chapter 8.5).

8.5 Unblocking the Cuff

Before unblocking the cuff, ensure that as little secretion as possible can enter the bronchi. While the cuff is being unblocked, suction secretions with a suction catheter passed through the tube.

If tubes with subglottic suction have been used, also suction the subglottic area before unblocking the cuff, see chapter 9.5.

To unblock the cuff, connect a syringe to the pilot balloon (**2a**). Completely remove all air/condensate if applicable. (See chapter 5. "General Precautions").

If it is not possible to unblock the cuff, apply the following measures:

- Remove the inner cannula. Make sure that the inflation line (**2a**) is not kinked and unblock the cuff again.
- If it is still not possible to unblock the cuff: carefully cut through the inflation line between the 15 mm connector and the neck flange (**5**).

8.6 Tube Change

Follow chapter 8.1 "Preparing the Tube".

.6.1 Tube Change for an Unstable Stoma

Insert a Seldinger wire into the trachea through the tracheostomy tube in situ. When removing the tube, approximately 3 hand's widths of the Seldinger wire should remain visible. The Seldinger wire remains in the trachea. Now insert the tube as per chapter 8.3.

8.6.2 Tube Change for a Stable Stoma

Assemble the white obturator (**15**) in a spare outer cannula. Then apply a pea-sized amount of lubricating gel to the tip of the obtu-

rator and tracheostomy tube to facilitate insertion.

Insert the obturator and the tube into the trachea through the tracheostoma. Then remove the obturator immediately, insert an inner cannula and block the cuff (as per chapter 8.3, item 4).

If necessary, a size 07 – 10 white obturator can also be used in combination with a Seldinger wire $\varnothing 1.27 \pm 0.04$ mm. The recannulation tube is suitable for all sizes.

9. Handling

9.1 Changing the Inner Cannulas

If viscous secretion collects in the inner cannula and cannot be suctioned, thus impeding airflow, replace the inner cannula with a new or cleaned inner cannula.

To remove the inner cannula, turn the 15 mm connector of the inner cannula anticlockwise (C).

Once a new inner cannula has been inserted into the outer cannula, turn the 15 mm connector clockwise until it locks into place, i.e. the markings are aligned with each other (C).

Inner cannulas with grooved low profile connectors function in the same way.

9.2 Fenestrated Tracheostomy Tubes

Unfenestrated inner cannulas have white 15 mm connectors and are used, for example, during mechanical ventilation.

To allow the patient to speak, insert a **fenestrated inner cannula** (blue 15 mm connector or blue grooved low profile connector) in the **fenestrated outer cannula** and unblock the cuff. Insert a speaking valve on the 15 mm connector of the inner cannula. Follow the instructions for use of the respective speaking valve.

The inner cannula REF 401 (not supplied) has an integrated speaking valve. The valve allows air to enter the trachea through the tube on inspiration; on expiration, the valve closes and the air flows alongside the tube, as well as through the fenestration, through the larynx into the upper respiratory tract, which allows vocalisation.

If there is an urge to cough or possible breathing difficulties, the REF 401 valve leaflets can be folded back and the patient can breathe unimpeded.

CAUTION: Speaking valves must only be used in patients who are awake and can breathe spontaneously. When a speaking

valve is introduced, the patient must be instructed by trained personnel and carefully monitored while becoming accustomed to it; make sure that the patient can breathe sufficiently. Unblock the cuff before using a speaking valve. Do not use a speaking valve during sleep.

CAUTION: Patients with laryngeal stenosis, vocal cord paralysis, severe tracheal stenosis, respiratory tract obstructions, respiratory tract infections or copious mucous secretions must not use a speaking valve.

9.3 Occlusion Cap of the Fenestrated Tube for Weaning

The fenestrated inner cannula can be closed by attaching the enclosed occlusion cap (10) to the 15 mm connector.

CAUTION: When preparing the tube, ensure that the patient's upper respiratory tract is unobstructed. The upper respiratory tract can be cleared by coughing out or suctioning the secretions. The cuff needs to be unblocked. When closing the tube, the patient's breathing and vital signs must be monitored by a nursing professional. If any signs of breathing difficulties are detected, the occlusion cap must be removed immediately.

9.4 Keeping the Fenestration Open

If the tube is used over a prolonged period, it must be ensured at regular intervals that the fenestration is not blocked by secretions, incrustation or in-growing tissue. If required, replace the tube.

9.5 Application of REF 306-P, REF 888-306-P with Subglottic Suction Line

These tubes are fitted with a flat suction line (9) that is fixed to the outer curve of the outer cannula and ends with two openings immediately above the cuff (2). The suction line has a female Luer connector, through which any secretion can be removed using a syringe. As an alternative, a special suction device with a vacuum controller connected to the enclosed connectors (11) can be used. After suctioning, ensure that the Luer connector is closed.

CAUTION:

- Ensure during suctioning that no excessive vacuum is used over a prolonged period (- maximum of 200 mbar).
- To prevent the subglottic area from drying out, we recommend the use of intermittent suction.
- The suction line may get blocked due to accumulated secretions or adhesions to the tracheal mucosa. If it is intended to

rinse the suction line (e.g. with a small amount of air or physiological saline solution), ensure beforehand that the cuff is sufficiently blocked (risk of aspiration).

- The suction line may cause pressure points in the area of the tracheostoma or may cause the tracheostoma to become irregular. In such cases, the doctor needs to decide whether this type of tube can still be used.

10. Cleaning, Disinfection and Storage

10.1 Cleaning

The tube (**1 + 3**) needs to be cleaned before reuse. Lukewarm potable water can be used for cleaning. TRACOE medical recommends the use of the cleaning products offered by TRACOE for cleaning. Once the tubes have been cleaned, rinse with potable water and air-dry. When cleaning the outer cannulas, take care not to damage the cuff.

CAUTION:

- Never use aggressive domestic cleaning agents, cleaning agents for dentures or solvents with a high alcohol content, for example, as these may have a lasting adverse effect on the function.
- Do not heat the tubes to a temperature above 65°C, as the safety of the products can then no longer be guaranteed.

10.2 Disinfection

Should disinfection appear necessary, TRACOE medical can provide a list of tested disinfectants on request. This list is also available at www.tracoe.com. The use of unsuitable disinfectants can result in damage to the product.

Once the tubes have been disinfected, rinse with sterile water and air-dry.

10.3 Storage

Store the clean tube in a dry location.

11. Product Modifications

Modifications to TRACOE products may only be carried out by employees or agents of TRACOE medical GmbH.

12. Disposal

These products may only be disposed of in accordance with the valid national regulations governing waste materials.

13. Returns

Returned used products will only be accepted subject to prior consultation and if a completed decontamination certificate is enclosed. This form is available either directly from TRACOE medical or via the website www.tracoe.com

14. General Terms and Conditions

The sale, delivery and return of all TRACOE products shall be effected exclusively on the basis of the valid General Terms and Conditions (GTC), which are available either from TRACOE medical GmbH or on our website at www.tracoe.com.

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